



DRILL – AIRE PLUS MINI INSTRUCTION MANUAL



Engler High Speed Dental Air Unit

COMPANY PROFILE

Engler Engineering Corporation has been in business since 1964 and occupies an 8000 square foot facility in Hialeah, Florida (USA). We manufacture ultrasonic dental scalers, polishers and combination units including electro surgery equipment and ultrasonic instruments for the veterinary market as well as a microprocessor controlled anesthesia delivery system and a respiratory monitor for veterinary use only.

We also manufacture dental equipment for the human market. Please visit our website www.englerusa.com for more detailed information or call us at the numbers shown below.

Engler Engineering Corp. acquired the exclusive manufacturing and marketing rights of **Dynax** products, including stretchers, animal restraint devices, comfort cots, heating pads, and other products. We also acquired the Alpha-Sonic, Ora-Sonic, and Pro-Sonic line of piezo scalers.

Engler Engineering Corporation's brand name veterinary products proudly include:

Excelsior, high speed dental air unit with vacuum / electro-surge / ultrasonic scaler / low speed / high speed / air / water syringe,

Son - Mate II, ultrasonic scaler / polisher,

Sonus II, ultrasonic dental scaler,

Poli - X, micromotor variable speed polisher,

Drill - Aire, high speed dental air unit, high speed, air / water syringe,

Drill - Aire Plus, high speed dental air unit, high speed, low speed, air / water syringe,

Scale - Aire Mini, high speed dental air unit with ultrasonic scaler / high speed / low speed / air / water syringe,

Scale - Aire, high speed dental air unit with ultrasonic scaler / high speed / low speed / air / water syringe and compressor / tank,

Tri - Mate, ultrasonic scaler / micromotor polisher / electro-surge,

A.D.S. 2000, microprocessor controlled anesthesia delivery system / ventilator,

Sentinel V.R.M., respiratory monitor.

Engler manufactures the **Sonus V** ultrasonic dental unit for the human market

We manufacture all of the inserts and tips used in the Engler products as well as many others on the market today in the 18K, 25K, and 30K frequency range.

Our repair department has the technical knowledge to repair and maintain most dental devices manufactured by other companies including Shorline.

Engler Engineering Corporation's foreign sales are handled through a large and growing network of dental and veterinary distributors. At the present time we are represented throughout Europe, South and Central America, Canada, Asia, New Zealand, Australia, the Middle East, and most other countries.

If you have any questions or comments, please contact:

Engler Engineering Corporation

1099 East 47th Street, Hialeah, Florida 33013

800-445-8581 – 305-688-8581 – FAX 305-685-7671

Web site: www.englerusa.com Help site: www.engler411.com

INTRODUCTION

Thank you for selecting the Drill - Aire Plus Mini High Speed Dental Air Unit. We believe you have selected the best product available for performing basic and advanced dentistry for your veterinary patients.

The design of the Drill - Aire Plus Mini uses time tested technology to produce a powerful and potent tool against periodontal disease.

The high speed drilling handpiece allows the operator to quickly and efficiently perform the same advanced dental techniques, drilling shaping and cutting to name a few, being taught in the largest teaching hospitals and clinics around the world. The low speed handpiece is used for smoothing and polishing the teeth after scaling.

PLEASE READ VERY CAREFULLY

Engler Engineering Corporation makes every effort to verify that all parts for the device along with any optional accessories ordered were shipped from our facility in Hialeah, Florida and are included in this shipment. It is imperative that you inspect the contents and if you find any pieces missing or damaged, you must notify us immediately. All claims submitted after fifteen days of receipt will not be valid.

All devices manufactured and / or sold by Engler Engineering Corporation are built and tested to approved standards. Any modification to the device, cables or hoses, initiated by others nullifies all warranty statements. Engler Engineering Corporation will not be held liable for any injury, death or damage of any type, due to non-authorized service and / or improper installation and / or improper use of the device.


This manual is not intended to teach dentistry. The information contained herein is intended only as a guide. Individuals not properly trained in dentistry should not use this equipment. It is intended for professional use only.

TABLE OF CONTENTS


COMPANY PROFILE.....	2
INTRODUCTION.....	3
INSTALLATION.....	5
DRILL - AIRE PLUS MINI DIAGRAM.....	6
BEFORE FIRST USE.....	7
HIGH SPEED HANDPIECE.....	7
LOW SPEED HANDPIECE.....	8
GETTING TO KNOW YOUR CONTROLS.....	10
WATER ON / OFF.....	10
AIR HANDPIECE TOGGLE SWITCH.....	10
HANDPIECE AIR FLOW CONTROL CONTROLS.....	10
WATER FLOW CONTROL.....	10
AIR HANDPIECE PRESSURE GAUGE.....	10
FOOTSWITCH CONTROL.....	10
Note: The high and low speed handpieces are NOT interchangeable.....	10
INSTRUCTIONS FOR OPERATION AND MAINTENANCE.....	11
REGULAR OPERATIONAL CONSIDERATIONS.....	11
MAINTENANCE KITS.....	11
START UP AND SHUT DOWN SEQUENCE.....	12
CLEANING HANDPIECES, FOOTSWITCH AND POWER CABLES.....	12
HIGH SPEED HANDPIECE.....	13
Cleaning the high speed handpiece.....	14
HIGH AND LOW SPEED HANDPIECE STERILIZATION PROCEDURES:	15
(AUTOClave AND CHEMIClave ONLY).....	15
IMPORTANT SAFETY PRECAUTIONS.....	16
INFECTIOn CONTROL.....	16
SGII – Doriot One - Piece Handpiece.....	18
THREE-WAY AIR / WATER SYRINGE.....	20
PROPHY ANGLE CLEANING AND MAINTENANCE INSTRUCTIONS.....	22
MAINTENANCE FREE PROPHY ANGLE.....	23
WATER TANK INSTRUCTIONS.....	24
WATER TANK CARE & MAINTENANCE.....	25
CARE AND MAINTENANCE OF YOUR PORTABLE WATER TANK.....	25
TECHNICAL SPECIFICATIONS:.....	26

INSTALLATION

1. Unpack unit and confirm all items checked off on your packing list..

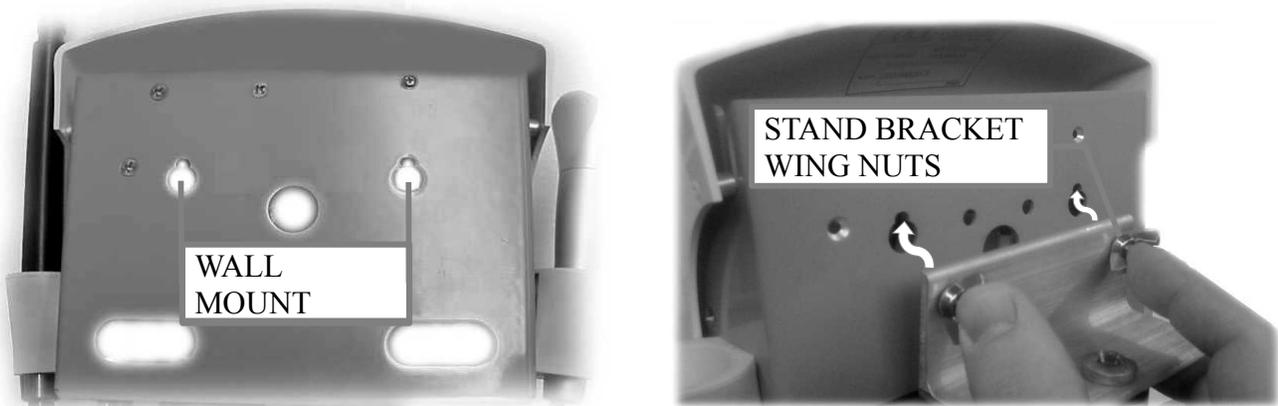

2. Attach unit to optional stand or wall mount (included):
 - A. Assemble stand.
 - B. Connect stand bracket to “key holes” on the back of the unit.
 - C. Tighten wing nuts.
3. Place handpieces on handpiece holders.

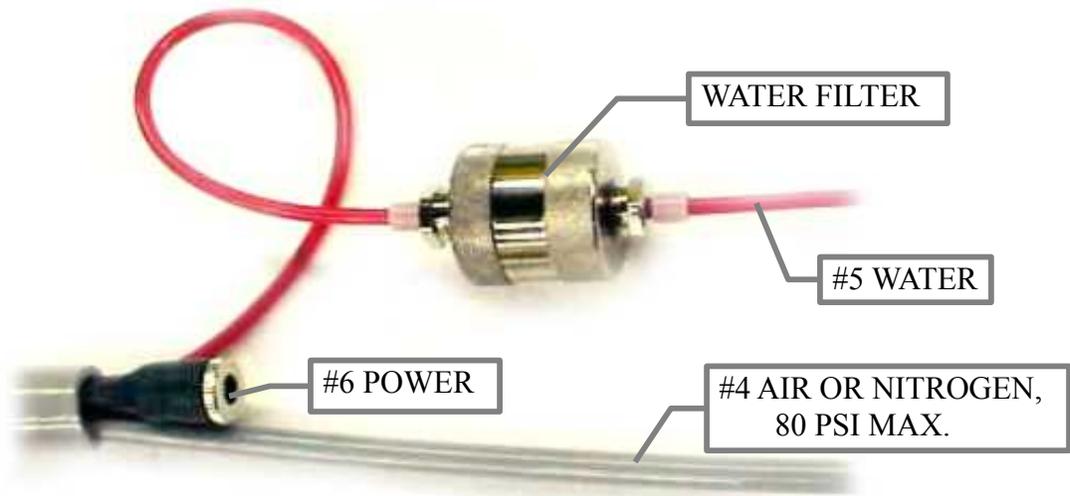


4. Connect pressurized air / nitrogen to gray hose. DO NOT EXCEED 80 PSI. Use pressure regulator and filter when needed. The recommended operating pressure is 75 PSI.
5. Connect water to red hose.
6. Connect power supply to power plug located at the end of the umbilical cord.
7. Connect power supply to power outlet.
8. Test each handpiece.

To avoid unit damage and / or injury, confirm input line pressure for the air or nitrogen is set to a maximum of 80 PSI. If using a compressor, verify your compressor pressure ratings, if it provides more than 80 psi, a pressure regulator is required.

DRILL - AIRE PLUS MINI DIAGRAM

Locate the handpieces on the front of your Drill - Aire Plus Mini. From right to left they are:


- High Speed Handpiece
- Low Speed Handpiece
- Air / Water Syringe

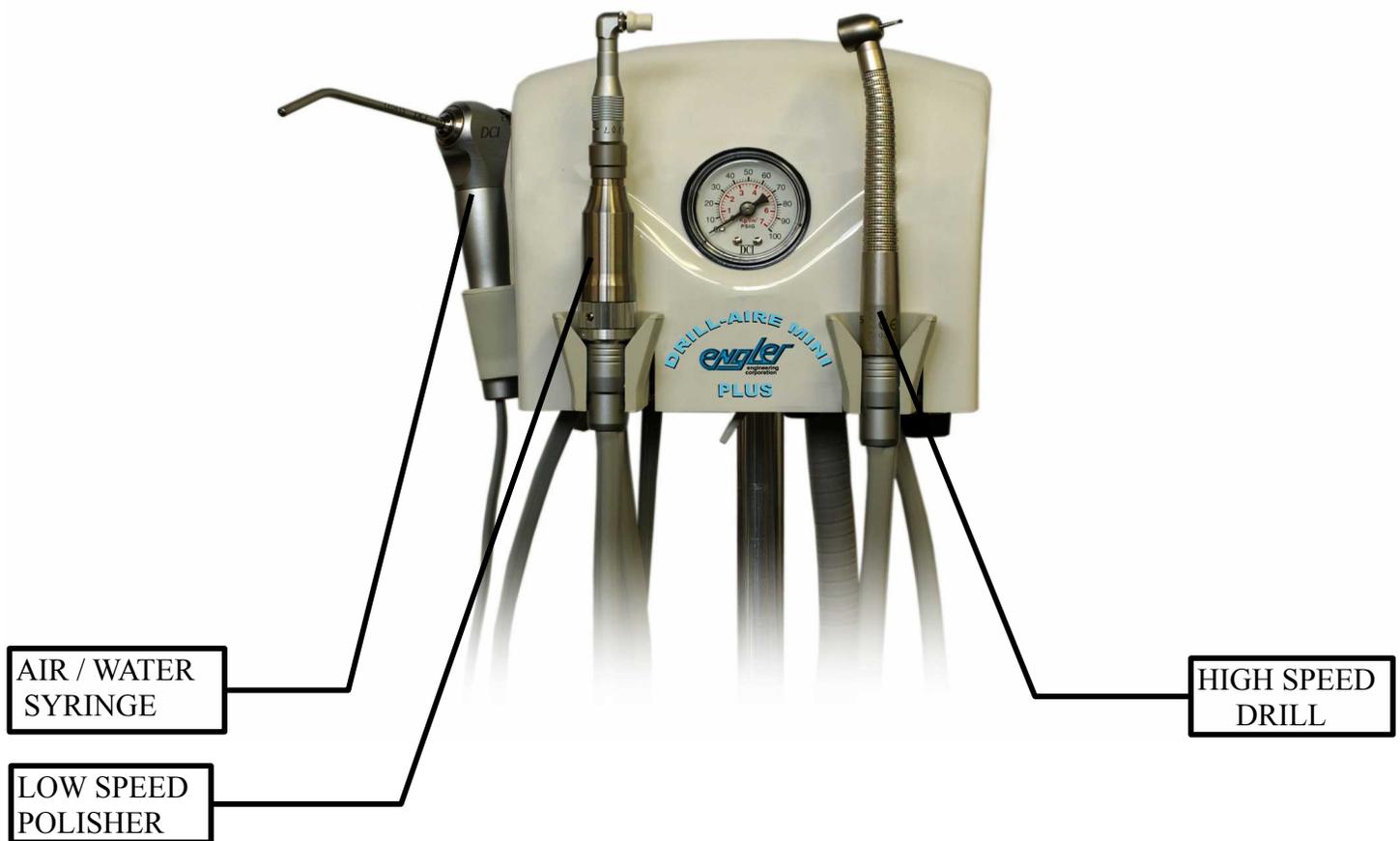

BEFORE FIRST USE

The Drill – Aire Plus Mini was purged before shipping, the first time the high speed handpieces are activated, the water lines will be empty. Turn the water regulator several rotations counterclockwise, lift the handpiece from it's cradle, holding the handpiece over a sink, turn the power ON by turning the power knob to click on, then press the footswitch. Water will flow after a few seconds.

HIGH SPEED HANDPIECE

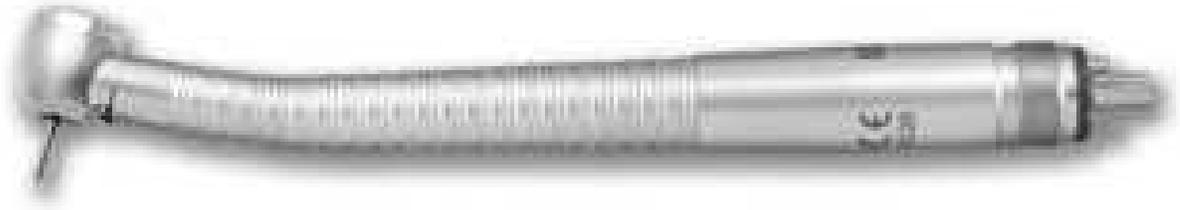


This handpiece is used for advanced dentistry, including but not limited to; cutting, sectioning, and shaping cracked or broken teeth, repairing, preparing cavities etc.



Flush to remove water. Keep well oiled.

LOW SPEED HANDPIECE

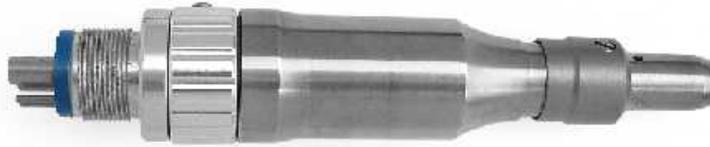

The Mini is equipped with the Doriot One-piece Handpiece.



Maximum rpm: 20,000

Attachments: Accepts both handpiece burs and Doriot / U-type attachments (prophy angle)
It is used with a prophy angle to polish the teeth after a scaling procedure. Use only approved attachments and polishing compounds. Follow all manufacturers recommendations.

The water should be turned OFF when using the low speed handpiece (polisher).

No straight handpiece is used with the short-grip Doriot one-piece handpiece since the prophy angle is attached directly to the motor.



Press-in here to release tip

AIR / WATER SYRINGE

The three-way air / water syringe features

- Well balanced design, and smooth styling for comfortable use.
- Easy release for a speedy exchange of tips.
- Fully autoclavable tips.

This handpiece allows the operator the ability to rinse the operative site with a stream of water or mist or dry / blow debris with a stream of air.



WATER ON / OFF

AIR FLOW CONTROL
(Behind Selector Switch)

LOW SPEED / HIGH SPEED HANDPIECE SELECTOR

FOOTSWITCH



GETTING TO KNOW YOUR CONTROLS

WATER ON / OFF

The control is equipped with a wet / dry toggle to activate the water flow. Move the toggle forward to turn water on.

AIR HANDPIECE TOGGLE SWITCH

This switch allows you to select between low speed and high speed handpieces when the scaler is not activated.

HANDPIECE AIR FLOW CONTROL CONTROLS

On the bottom side of the unit the flow control screws can be found. These screws control the amount of air that is delivered to the handpiece once the footswitch is pressed.

Caution: Allowing more pressure than the maximum allowed pressure to your handpiece will cause your handpiece to fail prematurely. Verify your handpiece manufacture's specifications before modifying these settings.

WATER FLOW CONTROL

This control adjusts the water flow to the handpieces. A water control knob is provided for each handpiece. Turn clockwise to decrease flow, and counter-clockwise to increase flow.

AIR HANDPIECE PRESSURE GAUGE

Gives a visual indication of the air pressure delivered to the handpieces.

FOOTSWITCH CONTROL

The foot control applies air pressure to the selected handpiece. The footswitch is responsive, high volume, and variable flow that can be actuated by pressing any point on the durable, chrome-plated cover. By pressing harder on the footswitch, more power will be delivered to the air handpieces.

Note: The high and low speed handpieces are **NOT** interchangeable

INSTRUCTIONS FOR OPERATION AND MAINTENANCE

REGULAR OPERATIONAL CONSIDERATIONS

The low speed polisher and the high speed drill must be lubricated daily when in use. Remove the quick disconnect hose and spray lubricant in the smaller of the two larger ports of the handpiece. The straight handpiece (if available) must be lubricated daily. Note: Short grip Doriot low speed handpiece does not use a straight handpiece.

The prophylaxis angle head must be removed, cleaned with soap and water, dried, and lubricated at the end of each day depending on use, then re-assembled.

All handpiece holders and three way syringe must be kept cleaned.

One can of Lares Handpiece Conditioner spray lubricant is included with every Drill – Aire Plus Mini.

For the best care and maintenance order the Drill - Aire Plus Mini deluxe maintenance kit which includes Lares Handpiece Conditioner from Engler Engineering Corporation (optional, not included with unit).

MAINTENANCE KITS

There are maintenance kits available for the Drill - Aire Plus Mini (optional, not included with unit). These kits are essential to keep your unit working properly. Each kit varies. Please make sure you specify your systems configuration and part number when ordering the kit.

KIT MINI BASIC , The standard kit, includes One Step conditioner for the low and high speed handpieces, and port cleaning tool necessary to perform basic maintenance.

KIT MINI DELUXE , includes replacement parts that are normal wear and tear items such as F.G. Bur kit, replacement O-rings for high and low speed handpieces in addition to items in the standard kit.

START UP AND SHUT DOWN SEQUENCE

Start up:

1. Connect unit to a clean, filtered water source.
2. Connect air / nitrogen source.
3. Set the air / nitrogen regulator to 75 PSI.

Shut down:

1. Clean all air handpieces according to the original manufacturer instructions. Thoroughly wipe the unit, all handpieces and cables with a mild cleaning solution or disinfectant and a damp cloth. Follow the procedures approved by your institution or use a validated infection control procedure. Do not allow fluids to enter the chassis. Do not autoclave the main unit.
2. Remove and sterilize three way syringe tip (tip only), prophy angle, burs, low and high speed handpieces. **Prophy angle, low and high speed handpieces must be lubricated regularly and after autoclave.**
3. Disconnect the air / nitrogen source.
4. Disconnect unit from the water source.

CLEANING HANDPIECES, FOOTSWITCH AND POWER CABLES

After each procedure, or at least once a day, it is suggested that the handpieces and cables be thoroughly cleaned and sterilized.

Clean the outer surface of the handpieces and cables with an antiseptic soap, rinse with water and sterilize with a chemical sterilization solution.

Caution: No chemicals or cleaners should ever be used inside or allowed to get into the scaler handpiece. Flush the handpiece thoroughly and completely with clean water.

The footswitch and power cables should be cleaned regularly by spraying a fine mist of sterilization or cleaning solution on the cables. It should remain on the cables for the length of time recommended by the manufacturer. Wipe the surface with a damp cloth and dry the cables completely.

HIGH SPEED HANDPIECE



INSERTING A BUR INTO THE HANDPIECE: CHANGING BURS

Lares high speed handpieces should be used with friction grip burs with shank diameters that conform to ISO and ADA size standards.

Model	Recommended	Optional
757 Ultralite / Euro	Standard (19.0 mm)	Surgical Length (26 mm) 757 Ultralite / Euro

1. Hold handpiece handle as shown and position end of thumb on push button with index finger wrapped around underside of handpiece neck for support.

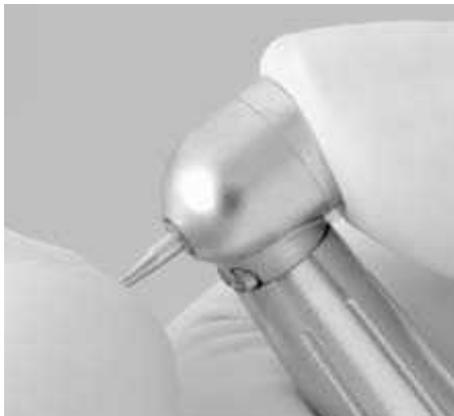


Figure 2

2. To insert a bur, first be sure bur is clean and free of external debris or corrosion. Without depressing push button, insert bur into chuck as far as possible. Then fully depress push button hard and fast while simultaneously inserting bur into chuck the rest of the way until fully seated. Release push button and insertion is completed.
Caution: Be sure to tug firmly on the bur immediately after completing the insertion procedure to verify full seating and secure retention of the bur before operation.
3. To remove a bur, fully depress push button hard and fast while simultaneously pulling bur until removed from chuck. **Figure 2.** Push button may then be released until next bur is inserted.

Caution: Never force a scored, bent or rusty bur into chuck as it may damage the chuck voiding the warranty. Never depress push button during handpiece operation or while turbine is still rotating. Be sure to remove bur at end of day.

Cleaning the high speed handpiece



1. Remove the bur.
2. Unscrew the high-speed handpiece from the air-hose coupler.
3. Use the water port cleaning tool (item # 10541) to clean the small water spray holes.
4. Use the spray lube (with red tube) to spray lubricant into the chuck and into the air drive hole (the smaller of the two large holes).
5. Re-connect the air-hose coupler.
6. Re-insert a bur or bur blank – NEVER run the handpiece without a bur or bur blank in place.
7. With water toggle switch OFF, depress the foot-switch pedal for at least 5 seconds to expel excess lube from the handpiece.
8. Wipe the handpiece of excess lubricant.

DAILY HANDPIECE LUBRICATION IS ABSOLUTELY ESSENTIAL.

Operate the high speed handpiece between 25 - 32 PSI. The handpiece is engineered to attain speeds of 360,000 RPM at 32 PSI.

Use the brush to remove foreign particles. A fine wire is provided for cleaning the water spray hole and to prevent clogging. Blowing air backward into the drill head can dislodge particles. NOTE: DO NOT attempt to blow particles from the rear end of the handpiece, as larger particles could block the water tube.

HIGH AND LOW SPEED HANDPIECE STERILIZATION PROCEDURES: (AUTOCLAVE AND CHEMICLAVE ONLY)

Clean external surface: Remove bur from handpiece and scrub with a toothbrush or 2" x 2" gauze using warm tap water. DO NOT IMMERSE HANDPIECE.

Dry: Thoroughly dry handpiece using gauze, paper towel or air syringe.

Clean / lubricate: Using a combination cleaner / lubricant, spray into handpiece drive air hole and in chuck.

Expel cleaner / lubricant: Reinsert bur into handpiece head, connect handpiece to tubing and run for a few seconds to thoroughly expel debris and excess lubricant.

Bag and cycle in autoclave: Place handpiece into autoclave bag or pouch. Cycle as per autoclave / chemiclave manufacturer's instructions. DO NOT EXCEED 275° F (135° C).

Cool down and lubricate: Allow handpiece to return to room temperature. Lubricate handpiece as per instructions listed above. Expel excess lubricant as per previous instructions. Wipe down handpiece.

DO

Use warm tap water to scrub the exterior of handpiece.

Expel excess lubricant from handpiece by running it for a few seconds after cleaning and lubricating.

Use separate cans of lubricant before and after sterilization to prevent cross contamination.

Use autoclave bags and pouches with indicators to protect handpiece.

DON'T

DON'T immerse handpiece in water or chemical disinfectant / sterilizer

DON'T use any type of disinfectant on handpiece.

DON'T sterilize handpiece with bur inserted.

DON'T exceed 275° F (135° C) in autoclave or chemiclave.

DON'T dry heat or heat transfer sterilize.

DON'T operate handpiece without bur or bur blank inserted in chuck.

TROUBLESHOOTING COMMON HIGH SPEED HANDPIECE ISSUES

A. Handpiece sluggish (could be lack of lubrication or too much debris)

1. Try heavily spraying "correct" drive air tube with a combination cleaner / lubricant #10083 available from Engler Engineering Corporation and run handpiece for a minute.

B. Burs are sticking inside of turbine or falling out

1. Flush the chuck / spindle with a handpiece cleaner.
2. Ensure that burs are not worn or scored, less than .0625" or larger than .0630" which can damage spindle.

C. Water spray is weak or completely stopped

1. Insert water port cleaning tool item # 10541 into water tube from head of handpiece to remove debris.
2. Use our Smart Cleaner to clear clogged tube.

If these simple solutions fail, more serious problems are likely affecting the handpiece, which should be sent to the manufacturer if under warranty or Engler Engineering Corporation if warranty has expired.

IMPORTANT SAFETY PRECAUTIONS

All high speed handpieces are potentially dangerous if safety precautions are not followed. Be sure to read and observe the following precautions.

Never use the back of the handpiece for tissue retraction, or otherwise cause push button to be depressed during operation. Doing so may result in button getting hot and burning the patient.

Never operate handpiece with a bent or damaged bur engaged in chuck.

Never insert or remove handpiece from oral cavity before rotation of bur is completely stopped.

Never operate handpiece at air pressure in excess of recommended maximum settings.

Never operate handpiece after turbine cartridge replacement without double checking that head cap is tightened securely.

Never depress push button during operation.

Never operate handpiece without fully inserting bur in chuck. Do not extend burs. Longer burs are available.

Sterilize between patients to prevent cross contamination.

INFECTION CONTROL

Be sure to use only Lares Handpiece Conditioner for this handpiece. Use of lubricants / conditioners other than Lares approved conditioner or failure to follow the maintenance schedule described above will automatically void the limited warranty for this product. Lares Handpiece Conditioner is available from Engler engineering corporation, part #10083.

All Lares high speed handpieces may be steam autoclaved or chemiclaved.

Prior to cycling, be sure to have Lares Handpiece Conditioner available with the appropriate nozzle hardware attached.

PROCEDURE

This procedure should be performed after every patient to prevent cross-contamination and to assure long, trouble-free operation.

1. Clean external surface of handpiece thoroughly to remove saliva, blood, and other organic soil. Scrub handpiece with small brush under running water. Rinse and dry thoroughly.

CAUTION:

NEVER Submerge Components In Any Cleaning Or Disinfecting Solution
DO NOT Use Ultrasonic Cleaners

2. Apply Lares Handpiece Conditioner. Follow specific instructions detailed on can.

3. Remove bur from chuck.

4. Place handpiece in autoclave bag. The use of autoclave bags prevents damage to the handpiece.

5. Load autoclave bag containing handpiece into autoclave or chemiclave. Be sure to load autoclave bags for maximum penetration of steam or chemical vapor.

6. Cycle the handpiece according to the autoclave / chemiclave manufacturer's instructions. Do not exceed 275°F (135°C).

CAUTION:

DO NOT Autoclave or Chemiclave For Extended / Unusual Periods Of Time (Such As Overnight).

DO NOT Leave Handpiece Components In Sterilizer After Cycle Is Completed.

7. Immediately remove handpiece from autoclave or chemiclave. Allow to cool sufficiently prior to handling.

SGII – Doriot One - Piece Handpiece

Your unit may be equipped with the **SGII – Doriot One-Piece Handpiece** instead of the **E-TYPE LOW SPEED HANDPIECE**, in this case use the following instructions.

The speed of the polishing head is proportional to the amount of pressure applied to the foot-switch pedal. Use low pressure to maintain a low speed.

SGII – Doriot One-piece Handpiece

Technical Facts


SGII: Doriot One-piece Handpiece

Maximum RPM: 20,000 min or 5,000 max

Attachment: accepts both handpiece burs and Doriot / U-type attachments

Air Requirements: Clean filtered moisture free air with recommended pressure of at least 45 psi. **Do not exceed 60 psi**.


Operation

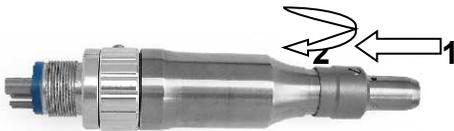

Removing handpiece bur or Doriot attachment from Doriot Handpiece: Hold the handpiece in the left hand, depress housing ring toward the body of the handpiece while twisting $\frac{1}{4}$ turn to the right to open the chuck.

****Push then Pull the handpiece bur or attachment and remove from the chuck.**

Forward / Reverse Speed:

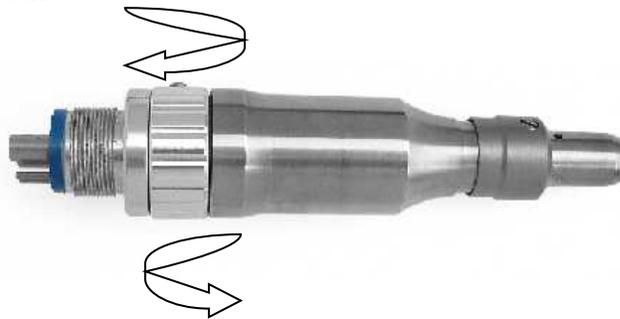

Forward:

Twist change ring fully clockwise

Reverse:

Twist change ring fully counter-clockwise

–In mid position, handpiece will not operate



Cleaning and Care

Your Doriot One-piece Handpiece low speed motor is a high quality precision instrument.

Incorrect maintenance and care will shorten the life of this product.

We recommend One Step Handpiece Cleaner and Conditioner for daily cleaning and lubrication, available from Engler Engineering Corporation.

External Cleaning

Remove angle or attachment from motor. Disconnect from hose and clean the external surface thoroughly with a sponge or gauze using warm tap water. **DO NOT IMMERSE INSTRUMENT.** Wipe dry with a clean cloth.

Clean / Lubricate

Lubricate after every sterilization. Using the combination cleaner / lubricant available from Engler Engineering Corporation. Spray the cleaner / lube into the drive air tube. Run motor to expell excess lubricant. This ensures all internal parts of motor are completed lubricated. Wipe off handpiece.



Weekly: Additional lubrication should be applied to inside nose of motor.

Sterilization:



Place motor into autoclave bag or pouch. Cycle per autoclave / chemiclave manufacturers' instructions. **DO NOT EXCEED 275° F (29 psi)/ 135°C (2 bar).**



CAUTION

DON'T use over-sized, bent, or scored burs.

DON'T immerse instrument in water or use ultrasonic cleaner.

DON'T use any type of disinfectant, chemical or soap on instrument. Use of chlorine products, aldehydes, etc. will damage handpiece and void all warranties.

DON'T operate motor while changing speed range, changing direction or inserting / removing attachments

DON'T exceed 275° F (29 psi) / 135°C (2 bar) in autoclave.

DON'T dry heat or heat transfer sterilize.

DON'T operate instrument without bur or bur blank inserted into chuck.



THREE-WAY AIR / WATER SYRINGE

This handpiece allows the operator to rinse the operative site with a stream of water or mist or dry / blow debris with a stream of air. The button on the left controls water flow. The button on the right controls airflow. Pressing both buttons at the same time provides mist. The air / water syringe works independently of the other handpieces and can be used alone or with the other handpieces.


The syringe features quick-change autoclavable tips: To remove a tip, press on the locking collar surrounding the tip socket and pull the used tip straight out of the socket. To insert a new tip, press locking collar and push tip into the socket as far as it will go. Release collar and gently tug on tip before using to ensure that tip is securely locked into socket.

SYRINGE TIP STERILIZATION

- 1) Remove contaminated syringe tip.
- 2) Remove all visible signs of contamination before autoclaving.
- 3) Autoclave tip at 132° C (270° F) for ten minutes.
- 4) Sterilize between patients.

NOTE: Since only the tips can be autoclaved, It is recommended that the air / water syringe be bagged with a disposable, single-use plastic sleeve between patients.

SEAL-TIGHT PROPHY ANGLE

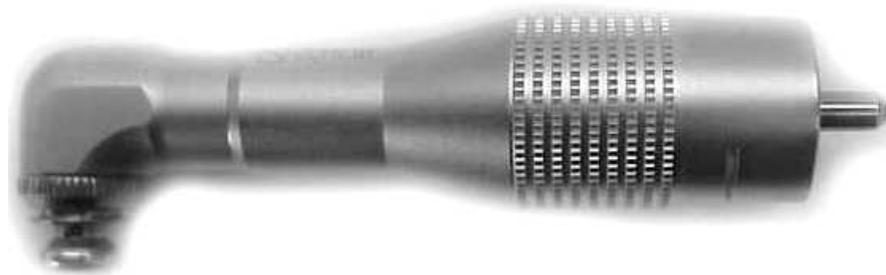

The Seal-Tight prophy angle comes standard with every unit.

This prophy angle is a precision-engineered dental device. All of the gears and shaft assemblies are made of high-grade stainless steel, which when kept clean and properly lubricated, will provide long, trouble-free service.

General Cautions;

When operating the prophy head always consider the safety of the patient.

Should the prophy head function abnormally, cease operation immediately. See below for maintenance instructions. If you need further assistance contact Engler Engineering Corporation for instructions.

Before use, always confirm that the brush or rubber cup is correctly seated in place.

If end cap unscrews by itself during a procedure; Switch the direction of the motor

Check that head and end cap are screwed together firmly. A poor fit between the head and end cap requires replacement of the entire prophy angle.

IMPORTANT: The prophy angle is rated for no more 3,000 R.P.M. – High speeds will result in the polished surface heating up, spattered polish, and shorter prophy angle life. Keep the unit set in the prophy range whenever using the prophy angle. **Always start with a low speed and then adjust to a higher speed as required.**

PROPHY ANGLE INSTRUCTIONS FOR USE

Dampen rubber cup.

Dab a small amount of polishing paste onto the rubber cup.

Place the prophy angle into the patients' mouth and lightly apply the rubber cup to the surface of the tooth with a circular motion.

Always begin the procedure using the lowest possible RPM, increase speed as necessary.

Do not allow the rubber cup to remain stationary on one area for an extended period of time. Keep the cup moving, so as to not overheat one area.

Add polish as needed. Move over all tooth surfaces and between teeth.

When finished, rinse the patient's mouth thoroughly with plenty of water.

PROPHY ANGLE CLEANING AND MAINTENANCE INSTRUCTIONS

The prophy angle is a precision engineered dental device. all gear and shaft assemblies are made of high grade stainless steel which must be kept free of debris. If cleaned and lubricated correctly will provide long, trouble-free service. The manufacturer recommends replacing prophy angles at least every 3 to 4 months depending on use. Prophyl angles may vary. Use the following instructions accordingly.

DAILY CLEANING AND LUBRICATION:




- 1.) Remove prophy angle from low speed handpiece.
- 2.) Discard used rubber cup.



- 3.) Remove head cap by turning **counterclockwise** to unscrew the knurled nut with the wrench provided.
- 4.) Wash the cap and head cavity thoroughly with a toothbrush in a bowl of warm soapy water.
- 5.) Rinse thoroughly with running water and shake off excess water.
- 6.) **DO NOT** attempt to dry this part with paper or cotton towels / swabs or gauze. Any particles left on the gears will keep them from turning properly. Use only alcohol to speed the drying process and / or a blow dryer to thoroughly dry the angle.
- 7.) Lubricate by placing one drop of mineral oil on the gears of the head cap and a drop inside the gear cavity.
- 8.) Being careful not to cross-thread, reassemble the prophy angle and wipe off all excess oil. Place a new rubber cup onto the head cap and confirm that the gears are meshing properly by rotating the cap – it should turn easily. If not, remove and try again. **DO NOT** use the wrench, only finger tighten.
- 9.) Slide the prophy angle onto the handpiece and lock the handpiece.

- 1.) Remove prophy angle from low speed handpiece.
- 2.) Discard used rubber cup.

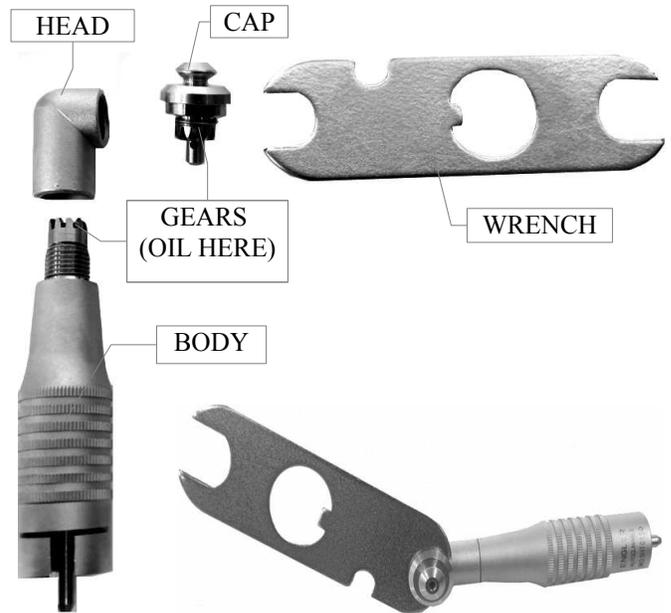

- 3.) Use the wrench to remove the cap from the head.
- 4.) Turn **clockwise** to unscrew the head (top portion) from the body (bottom portion).
- 5.) Place the cap, head and body into a bowl of hot soapy water.
- 6.) Wash thoroughly with a toothbrush.
- 7.) Rinse well with running water and shake off. **DO NOT** attempt to dry these parts with paper or cotton towels / swabs or gauze. Any particles left on the gears will keep them from turning properly. Use only alcohol to speed the drying process and / or a blow drier to thoroughly dry the angle.
- 8.) Lubricate by placing one drop of mineral oil on each gear (see diagram).
- 9.) Being careful not to cross-thread, reassemble the prophy angle and wipe off all excess oil. Place a new rubber cup on the onto the cap and confirm that the gears are meshing properly by rotating the cap – it should turn easily. If not, remove the cap and try again.
- 10.) Slide the prophy angle onto the handpiece and lock the handpiece.

OPTIONAL ACCESORIES

MAINTENANCE FREE PROPHY ANGLE



P-MF

P-MF Maintenance free prophy angle



P-106

P-106 screw on rubber cups 144/pkt

CARE and STERILIZATION PROCEDURES

After each prophy:

Rinse abrasive paste from head and cup area with water. Then remove cup.

Thoroughly clean the outside of angle with disinfectant.

Place angle into a sterilization bag

Follow sterilizer manufacturer's recommendations.

Do not exceed 275 ° F(135 °C).

Keep angle in bag until ready for use.

You are now ready for your next prophy.

CAUTIONS AND WARNINGS:

Sterilize prior to disposal

Do not attempt to disassemble.

DO NOT submerge in liquids, including ultrasonic solutions.

Operate handpiece in the forward direction (counterclockwise when facing you) to prevent threaded prophy cups from unscrewing during the procedure.

Not recommended for use above 3000 R.P.M.

If the head of the angle becomes hot during use, lubricate the rim of the cup / screw hole with mineral oil.

Do not heat over 275 °F (135 °C).

Use only Engler Care Free Prophy Rubber Cups. Other brands will not properly seal the angle, causing premature wear and voiding the warranty.

Use 1 year, sterilize, then dispose of properly.

YOUR CAREFREE ANGLE IS WARRANTED AGAINST DEFECTS IN MATERIAL AND WORKMANSHIP FOR 6 MONTHS. A COPY OF OUR INVOICE OR PICKING TICKET WILL BE REQUIRED AS PROOF OF PURCHASE.

CAREFREE ANGLE WARRANTY IS VOID IF

Engler Care-Free rubber cups are not used exclusively.

Sterilization procedure is not followed properly.

The angle has been submerged in any liquid.

The angle has been damaged or abused.

Damaged due to use at high speed.

WATER TANK INSTRUCTIONS



PT-1 Water Tank

DIRECTIONS:

1. Remove pump and cap assembly.
2. Fill tank with distilled water or medicated solution up to the "FILL LINE". Do NOT fill beyond this line.
3. Replace pump and cap assembly and tighten securely.
4. Pressurize tank by pumping it approximately 20 - 40 times (depending on the amount of liquid used). If a hissing sound is detected, tank is over-pressurized. Stop pumping. Leave tank on a level surface until hissing stops. Insert the male quick disconnect on the end of the waterline from the Drill - Aire Plus Mini into female quick disconnect provided on tank.
5. Release air pressure by PULLING and TURNING pressure relief valve, located on the side of the bottle.

WATER TANK CARE & MAINTENANCE

- A. Release air pressure by pulling and turning knob of pressure relief valve. Pull out fully and allow air to escape.
- B. Loosen cap slowly. Remove pump & cap assembly. Pour out any remaining liquid & rinse all parts thoroughly with clean water.
- C. Always store tank empty and with tank cap loose.

TROUBLESHOOTING:

PROBLEM: TANK FAILS TO PRESSURIZE.

- Be sure cap is tight.
- Check to see if pressure relief valve is closed.
- Remove the pump from the tank. Turn pump handle counterclockwise and lift handle to unlock. At the top of the pump cap there is an opening that indicates "oil here". Place 3 - 5 drops of mineral oil into the opening. Pump several times to work the oil into the walls of the pump until it moves freely. Repeat if necessary. Screw the pump back into the tank and resume normal operations. This process should be repeated often as necessary and depending on usage, or when pump starts to work harder.
- Black particles found in water bottle indicates that the pump assembly is deteriorating. Order new pump assembly from Engler Engineering.Corp.

Pump assembly has been pre-lubricated prior to shipping.

WARNING:

READ AND FOLLOW ALL INSTRUCTIONS.

ALWAYS INSPECT your pump before each use.

ALWAYS RELEASE AIR pressure before removing pump or servicing tank, by pulling pressure relief valve knob out fully.

DO NOT use mechanical devices to pressurize the tank .They can create excessive and dangerous pressure which could cause the tank to explode.

DO NOT STAND over pressurized tank while using it or pumping it

DO NOT USE solutions warmer than 105F.

DO NOT damage or alter the functions of the pressure relief valve or plug the pressure relief valve hole, as this could cause the tank to explode

DO NOT pressurize the tank until ready for use.

DO NOT lift or carry the tank by waterline, extension rod or pump handle unless it is securely locked in place.

CARE AND MAINTENANCE OF YOUR PORTABLE WATER TANK

TO KEEP SLIME FROM FORMING INSIDE THE TANK AND EVENTUALLY GETTING INTO THE DENTAL UNIT CAUSING IRREVERSIBLE DAMAGE:

1. Every two weeks dispose of water in tank. Pour ½ gallon of hot water and 1 ounce bleach into tank and swirl the liquid thoroughly inside the tank.
2. **Dispose of bleach mixture and rinse tank with clean water thoroughly and completely.**
3. Clean the outside of the pump / tank according to your facilities normal cleaning procedures.
4. The pump assembly has been pre-lubricated. DO NOT TAKE THIS ASSEMBLY APART.

WARNING: TAKE CARE TO RINSE TANK THOROUGHLY AND COMPLETELY AFTER CLEANING. ANY CLEANERS REMAINING IN THE TANK COULD DAMAGE INTERNAL COMPONENTS OF THE DRILL - AIRE PLUS MINI AND VOID THE WARRANTY.

TECHNICAL SPECIFICATIONS:

Input pressure (air or nitrogen):

Operating 75 psi
Maximum 80 psi
Minimum 40 psi

Low speed handpiece:

20,000 RPM (5000 RPM for prophy)

High speed handpiece:

360,000 RPM

Power supply:

Input: 100-240 V ~, 1.5A, 50 - 60Hz
Output: 24 VDC, 2.5 Amps

Control box dimensions:

9" X 4" X 6"

Net weight:

8 lbs.

Shipping weight:

14 lbs.

Shipping box dimensions:

14" X 14" X 14"

Engler Engineering Corporation
1099 East 47th Street, Hialeah, Florida 33013
800-445-8581 – 305-688-8581 – FAX 305-685-7671
Web site: www.englerusa.com Help site: www.engler411.com